## High Performance Cooled CCD Camera System ALTA U6




The Alta U6 uses a large format megapixel Kodak Blue Plus sensor. This CCD has two output amplifiers: one for low noise and a second for exceptionally high dynamic range. Please specify which amplifier is desired when placing an order.

Imaging Area of CCD


- 1024 x 1024 array, 24 x 24 micron pixels
- 5 MHz 12-bit and 1 MHz 16-bit digitization
- 32Mbyte camera memory
- USB 2.0 interface: no plug in cards or external controllers
- Programmable, intelligent cooling to 50°C below ambient
- Binning up to 8 Horizontal x 1024 Vertical
- Subarray readout and fast sequencing modes
- Precision time delayed integration (TDI) and kinetics mode readout
- Programmable fan speed for low / zero vibration
- Two serial port outputs for control of peripheral devices
- General purpose programmable I/O port
- External triggering and strobe controls
- ActiveX drivers included with every system
- Field upgradeable firmware
- Fused silica windows
- Runs from single 12V supply with input voltage monitor
- Compact enclosure
- Programmable status indicators



- Radiology
- Optical testing
- Non-destructive testing

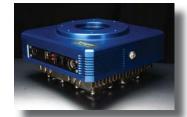



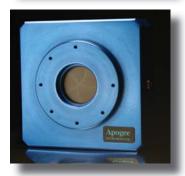



## **CCD** SPECIFICATIONS

CCD	Kodak KAF-1001E
Array Size (pixels)	1024 x 1024
Pixel Size	24 x 24 microns
Imaging Area	24.6 x 24.6 mm (604 mm2)
Imaging Diagonal	34.8 mm
Video Imager Size	2.17"
Linear Full Well (typical)	Low noise output: 200K electrons
Ellicar Full Well (typical)	High DR output: 500K electrons
Dynamic Range	<u> </u>
	High DR output: 500K electrons
Dynamic Range	High DR output: 500K electrons >87 dB

For complete CCD specifications, including cosmetic grading, see data sheet from manufacturer.






ALTA U6 era System Performance	
	ZILIIZ

PC Interface	USB 2.0
Max. Cable Length	5 meters between hubs; 5 hubs maximum (max. total of 30m)
Digital Resolution	16 bits at 1 MHz and 12 bits at 5 MHz
System Noise (typical)	8 e <sup>-</sup> RMS at 1 MHz and 2 counts at 5 MHz
Pixel Binning	1x1 to 8 x 1024 on-chip
Exposure Time	30 milliseconds to 183 minutes (2.56 microsecond increments)
Image Sequencing	1 to 65535 image sequences under software control
Frame Sizes	Full frame, subframe, focus mode
Cooling (typical)	Thermoelectric cooler with forced air. Maximum cooling 50°C below ambient temperature
Dark Current (typical)	1 e <sup>-</sup> /pixel/sec (-25°C)
Temperature Stability	$\pm 0.1$ °C
Camera Head Size	D2. Low profile: D6. Aluminum, hard blue anodized. 6" x 6" x 2.5" (15 x 15 x 6.35 cm) Weight: 3.1 lb. (1.4 kg)
Mounting	3.5" bolt circle. 2" 24 tpi thread. Optional Nikon F-mount or Canon FD mount.
Back Focal Distance	Standard: 1.020" (25.91 mm). Low profile: 0.455" (11.56 mm). [optical]
Operating Environment	-22° to 27°C. Relative humidity: 10 to 90% non-condensing.
Cable Length	Standard: 15 ft (4.5m)
Power	40W maximum power with shutter open and cooling maximum. AC/DC "brick" supply with int'l AC input plug (100-240V, 50-60 Hz). Alternate 12V input from user's source.
Shutter	Standard: Melles Griot 43mm. Low profile: no shutter.
Remote Triggering	LVTTL input allows exposure to start within 25 microseconds of rising edge of trigger

## **CCD SENSITIVITY**

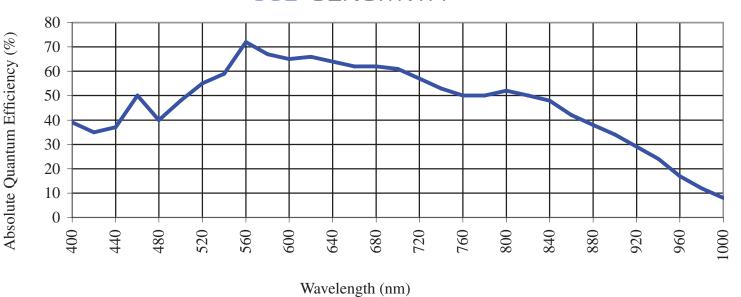

©2007 Apogee Instruments Inc. Alta is a registered trademark of Apogee Instruments Inc. Specifications subject to change without notice.

